

# Seeing is Believing: Innovative Use of Raman Microscopy to Visualize the Integrity of a Barrier Membrane Coating for Taste-Masking of Acetaminophen Granules in Chewable Tablets

Raxit Y. Mehta, Charles R. Cunningham, David Ferrizzi and Ali R. Rajabi-Siahboomi

Colorcon, Inc., Harleysville, PA USA

AAPS  
Poster Reprint 2016

## Purpose

Effective taste-masking of particles for incorporation into chewable dosage forms is influenced by the consistency and integrity of the applied barrier membrane (BM) coating. During the tablet compression process, the coating around drug particles may fracture leading to an adverse taste profile from the final dosage form. In addition to the barrier coating composition and application levels, the substrate morphology, film strength, and particle size also influence the final integrity of the coated particles.<sup>1,2</sup> In this study, Raman microscopy was used to generate images and visualize the integrity of barrier membrane coatings on two different grades of acetaminophen (APAP) granules with differing morphologies. Raman imaging was used to evaluate granule characteristics from two stages of the manufacturing process; after coating and after compression of the coated granules. Further, the observations from the Raman image analysis were correlated with in vitro drug release.

## Experimental Methods

### Coating of APAP Granules

For taste-masking purposes, two grades of APAP, Compap (amorphous, spray dried, d(50): 181  $\mu$ m) and Special granular (crystalline d(50): 332  $\mu$ m) (Covidien, USA) were coated using Surelease<sup>®</sup> Ethylcellulose Dispersion Type B NF, combined with a hypromellose-based Opadry<sup>®</sup> complete film coating system (Colorcon Inc., West Point, USA) used as a pore-former at a ratio of 85:15 w/w. The coating was applied to the granules using a top spray fluid bed coater (Glatt GPCG-2, USA) up to 30% weight gain (WG).

### Chewable Tablet Formulation

The coated APAP granules were passed through a 20 mesh sieve, then blended with the Parateck ODT blend (Merck Millipore, DE) plus a sweetener, disintegrant and colloidal silica. Magnesium stearate and FD&C blue #1 aluminum lake were passed through a 60 mesh sieve, then added to the chewable tablet blend as lubricant and colorant. The chewable tablet blend was compressed using a single station manual compression press (Globe Pharma, USA) with 12.5 mm flat-faced beveled edge tools at a compression pressure of 1200 psi. The target tablet weight was kept constant at 770 mg. The composition of the chewable tablet formulation is shown in Table 1.

**Table 1. APAP Chewable Tablet Formulation**

Ingredients	Supplier	(%)
Coated APAP granules	--	13.9
Parateck ODT blend	EMD Millipore, USA	76.8
NutraSweet (Aspartame)	NutraSweet, USA	0.8
Kollidon <sup>®</sup> CL-F (Crospovidone)	BASF, DE	5.0
Cab-o-Sil M5-P (Colloidal silica)	Cabot Corp., USA	1.5
Magnesium stearate	Peter Greven GmbH, DE	1.8
FD&C Blue #1	Colorcon, USA	0.2
<b>Total</b>		<b>100.0</b>

## Raman Spectroscopy Studies

The Raman maps of coated APAP granules and chewable tablets were acquired using a Horiba LabRam HR at laser excitation wavelength of 785 nm and power <40 mW. For Raman microscopic analysis, samples of the coated APAP granules were entrapped in UV cure epoxy resin and the cross-section obtained using microtome knife; whilst cross-section for the chewable tablets was directly obtained using the microtome knife.

## Dissolution Studies


In vitro dissolution testing of the coated granules and chewable tablets was carried out using USP Apparatus II (paddles) at 75 rpm in 900 ml of pH5.8 phosphate buffer. Drug release was determined spectrophotometrically at a wavelength of 243 nm.

## Results

### Raman Microscopy of Coated Granules

Raman images of the coated granules confirmed the presence of the coating surrounding the particle surfaces for both grades of APAP (Figure 1). However, non-uniformity of the coating around the edges of the irregularly shaped Compap granules was clearly observed. Compared to the Compap granules, the coating on the more regular shaped Special granular grade appeared more uniform in thickness, continuity, and integrity.


**Figure 1. Raman Imaging of Coated APAP Granules prior to Compression**

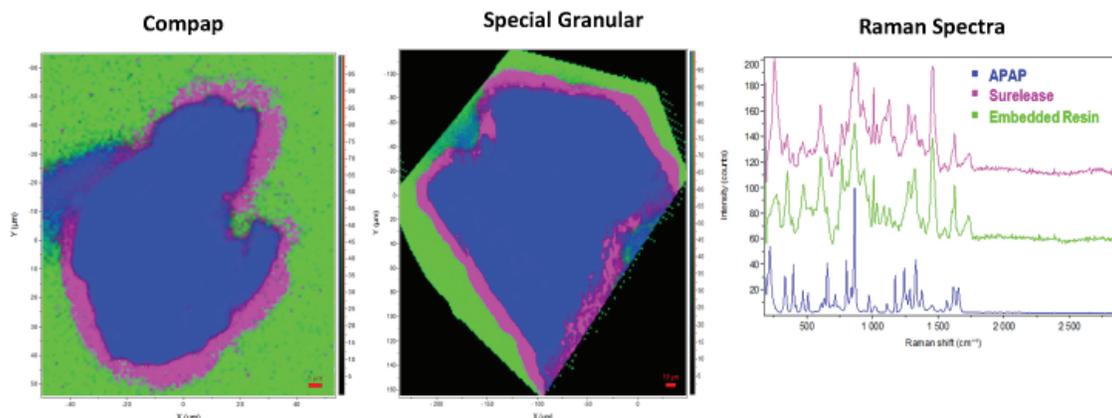

### Raman Microscopy of Compressed Chewable Tablets


The Raman spectra of the compressed tablets indicate clear differentiation of the individual formulation components (Figure 2).

**Figure 2. Raman Spectra of APAP Compressed Chewable Tablets**

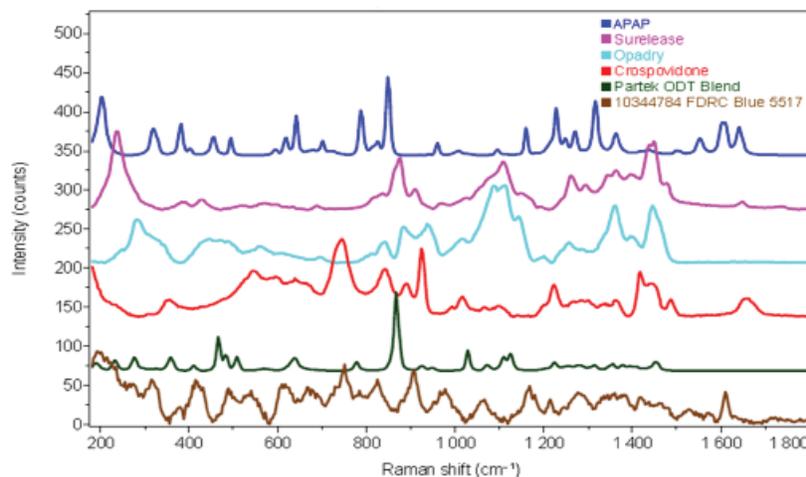


In chewable tablets, the integrity of BM coating layer around the Compap particle was further disrupted, probably due to the effect of compression force. The coating layer (highlighted in purple) was diffused into the other regions of the tablet leaving the Compap particle unprotected. Conversely, the Special granular particles coating layer was found to be intact and adhered uniformly (Figure 3).

**Figure 3. Raman Imaging of Single APAP Particle after Compression**

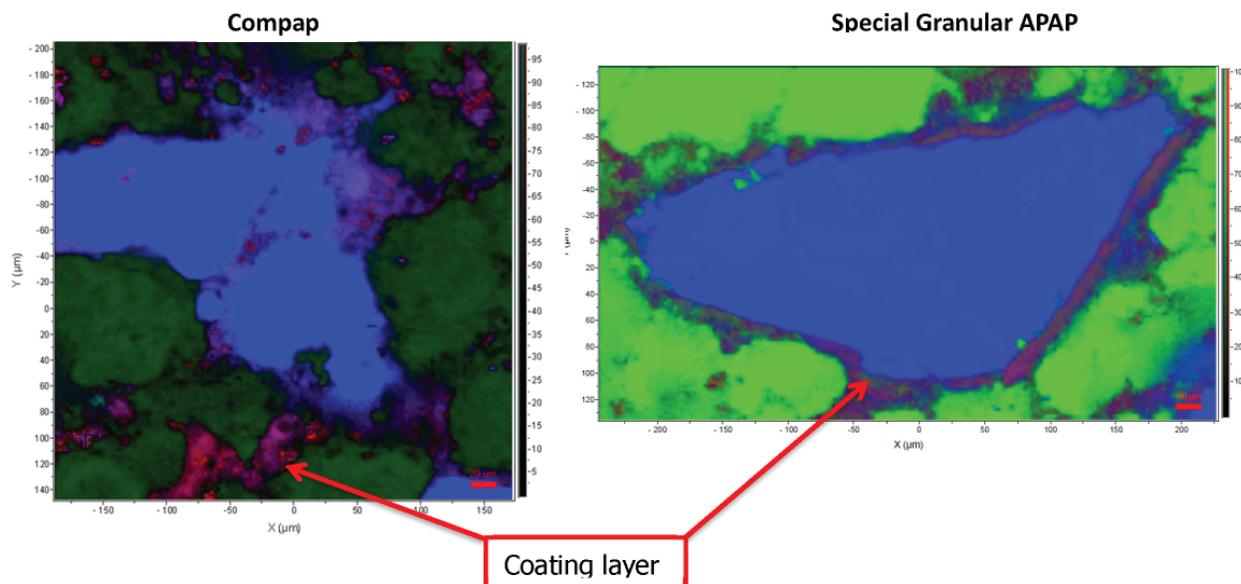

Further individual component spectral analysis was carried out by retaining the coating layer (Surelease) component spectra (purple) and eliminating all other component spectra from the analysis. This highlighted the differences in coating layer integrity between the particles of different morphology (Figure 4).


**Figure 4. Raman Imaging of Surelease Coating around APAP Particle within Chewable Tablet**

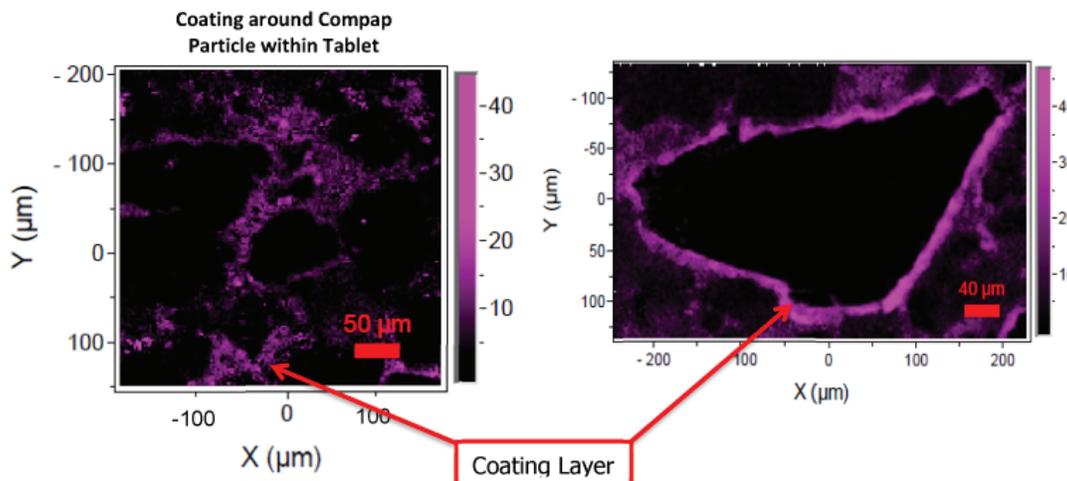


### Drug Release Analysis

The observed differences in coating layer integrity after compression of the coated granules was further illustrated by dissolution testing. Slower initial drug release, in the first 5 mins, is necessary to attain efficient taste-masking of acetaminophen granules. The drug release from coated Compap granules was significantly faster than Special granular grades, primarily due to non-uniformity of coating around the irregular particle shape. Further, the compression of the Compap granules into chewable tablets resulted in a 41% increase in drug release in the first 5 mins. Whereas, compression of Special granular grades resulted in only a 5% rise in drug release after 5 mins of dissolution (Figures 5 and 6).

**Figure 5. Coated Granule Drug Release Before and After Compression**

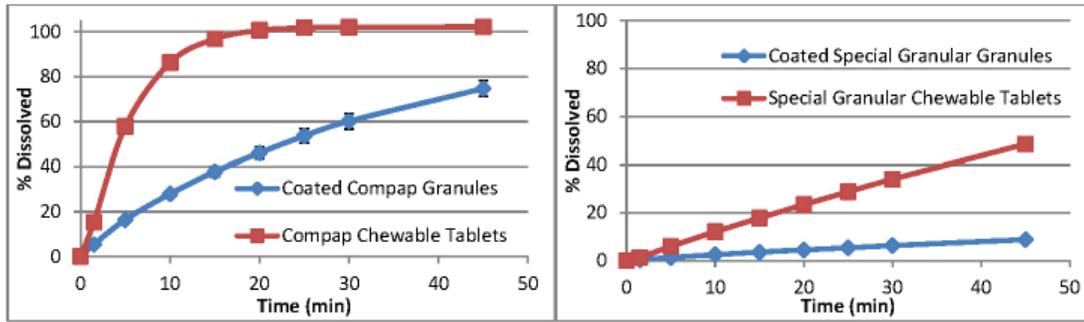

**Figure 6. Expanded View of Coated Granule Drug Release Before and After Compression**

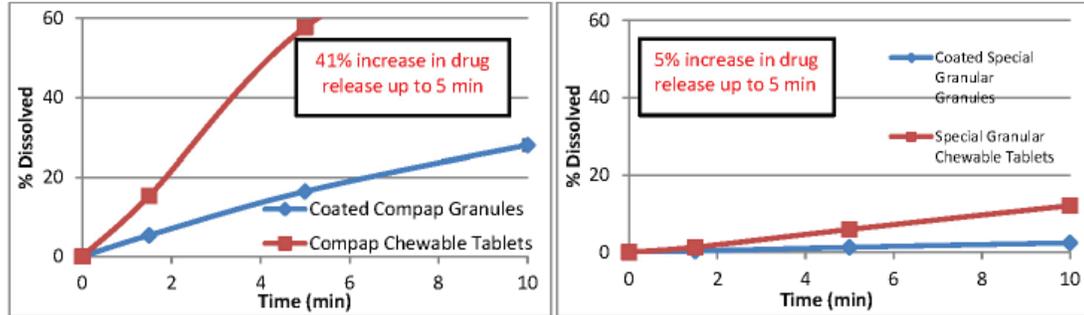

Once substrate considerations are understood, the drug release profile can be tailored to achieve either greater lag times or faster release by varying the coating weight gain and Surelease:Opadry ratios.<sup>3</sup>

## Conclusions

In this study, drug particles were successfully coated with Surelease:Opadry giving a reduced initial drug release that may provide desirable taste-masking characteristics when compressed into chewable tablets. The effect of substrate morphology and particle characteristics were found to be important in maintaining the integrity of the coating throughout the process of manufacturing chewable tablets. Raman microscopy was a useful technique to visualize the integrity of the coating layer at various stages of chewable tablets development and provided evidence for APAP release profile patterns seen from their compressed chewable tablets.

## References

1. Pearnchob N, Bodmeier R., Coating of pellets with micronized ethylcellulose particles by a dry powder coating technique. *International Journal of Pharmaceutics*. 2003; 268:1-11.
2. Terebesi I, Bodmeier R., Optimised process and formulation conditions for extended release dry polymer powder-coated pellets. *European Journal of Pharmaceutics and Bio pharmaceutics*. 2010; 75:63-70.
3. Vesey C, Farrell T, Rajabi-Siahboomi AR., Evaluation of alternative plasticizers for Surelease®, an aqueous ethylcellulose dispersion for modified release film-coating. *Controlled Release Society Annual Meeting* (2005).

## Acknowledgements

The authors are grateful to Mr. Mark Sparrow of the RJ Lee Group for their support of the Raman Microscopy testing and analysis.

The information contained herein, to the best of Colorcon, Inc.'s knowledge is true and accurate. Any recommendations or suggestions of Colorcon, Inc. with regard to the products provided by Colorcon, Inc. are made without warranty, either implied or expressed, because of the variations in methods, conditions and equipment which may be used in commercially processing the products, and no such warranties are made for the suitability of the products for any applications that you may have disclosed. Colorcon, Inc. shall not be liable for loss of profit or for incidental, special or consequential loss or damages.

Colorcon, Inc. makes no warranty, either expressed or implied, that the use of the products provided by Colorcon, Inc., will not infringe any trademark, trade name, copyright, patent or other rights held by any third person or entity when used in the customer's application.

For more information, contact your Colorcon representative or call:

North America  
+1-215-699-7733

Europe/Middle East/Africa  
+44-(0)-1322-293000

Latin America  
+54-1-5556-7700

India  
+91-832-672373

China  
+86-21-61982300

Surelease®

You can also visit our website at [www.colorcon.com](http://www.colorcon.com)



BPSI Holdings, LLC 2016

All trademarks, except where noted, are property of BPSI Holdings LLC. The information contained in this document is proprietary to Colorcon, Inc.

pr\_aaps-sur\_raman\_chewable\_11\_2016