

Evaluation of Formulated Ethylcellulose Aqueous Dispersions to Develop Controlled Release Multiparticulate Formulation of Freely Soluble Drug, Propranolol Hydrochloride

Kevin Hoxha, Raxit Mehta, Manish Ghimire, Ali Rajabi-Siahboomi
Colorcon, Inc. 275 Ruth Road, Harleysville, PA, USA
www.colorcon.com

CRS
Poster Reprint 2021

Purpose

In this study, the controlled release multiparticulate formulation of a freely soluble model drug, propranolol HCl, was evaluated using a fully formulated coating Surelease[®], Ethylcellulose Dispersion and Aquacoat[®] ECD, an ethylcellulose dispersion that allows for the addition of desired type and quantity of a plasticizer.

Methods

Propranolol hydrochloride was layered onto sugar spheres (850 – 1000 µm) (Suglets[®], PF011 Colorcon, USA) using an HPMC-based Opadry[®] (Colorcon, USA) as binder with 70:30 w/w drug/binder ratio in Bosch Hüttlin Unilab. The multiparticulates (MP) were then barrier membrane coated either with Surelease (Colorcon, USA) or with Aquacoat ECD-30 (Colorcon, USA) in a fluid bed coater equipped with Wurster configuration (Glatt GPCG-2). Before coating, a water soluble plasticizer, triethyl citrate (TEC, Vertellus, USA) at 24 % w/w with respect to polymer, was added to Aquacoat ECD. Ethylcellulose aqueous dispersions were coated to a 10 % w/w theoretical weight gain (WG). A top-coat (using the same Opadry) at a 1% w/w WG was applied over the barrier membrane coatings. Each batch of coated MPs was split into three equal parts, then two parts were subjected to post coating curing for 2 hours and 24 h in a 60°C oven. The drug release from MPs was characterized at 289 nm wavelength using USP Apparatus II (paddles), 1000 mL DI water at 37°C. The dissolution results generated were compared using the f_2 factor.¹ An f_2 value between 50 and 100 indicates that the two dissolution profiles are similar. The morphology of the barrier membranes was examined using SEM (Phenom XL, Thermo Fisher Scientific, USA).

Table 1: Batch Preparation Overview

Formulation Stages	Component	Solids (%w/w)	Weight Gain (%w/w)
Drug layering	Propranolol HCl Opadry	15	7
Barrier membrane coating	Surelease	15	10
	Aquacoat ECD-30 24% w/w TEC	15	10
Top-coat	Opadry	10	1

Table 2: Suglets Coating Parameters

Process Parameter	Drug Layering	Ethylcellulose Coating
Equipment	Huttlin Fluid Bed	Glatt CPCG-2
Batch size (kg)	4 – 6 kg	2.0
Inlet temperature (°C)	57 - 58	53 - 65
Product temperature (°C)	45 - 50	41 - 46
Outlet temperature (°C)	41 - 46	38 - 44
Atomizing air (bar)	1.5	2.0
Air volume (CMH) / (m ³ /hr)	325	120 - 130
Nozzle size (mm)	1.2	1.2
Fluid delivery rate (g/min)	12 - 24	10 – 14

Results: Dissolution

Stable drug release performance, independent of curing duration, was observed for the Surelease coated multiparticulate systems.

Aquacoat ECD coated multiparticulates required a minimum of two hours of curing to ensure full coalescence of the film and consistent slow drug release.

Figure 1. Drug Release from Propranolol HCl MR MP, Coated with Fully Formulated Ethylcellulose Dispersion (Surelease); Effect of Curing Duration



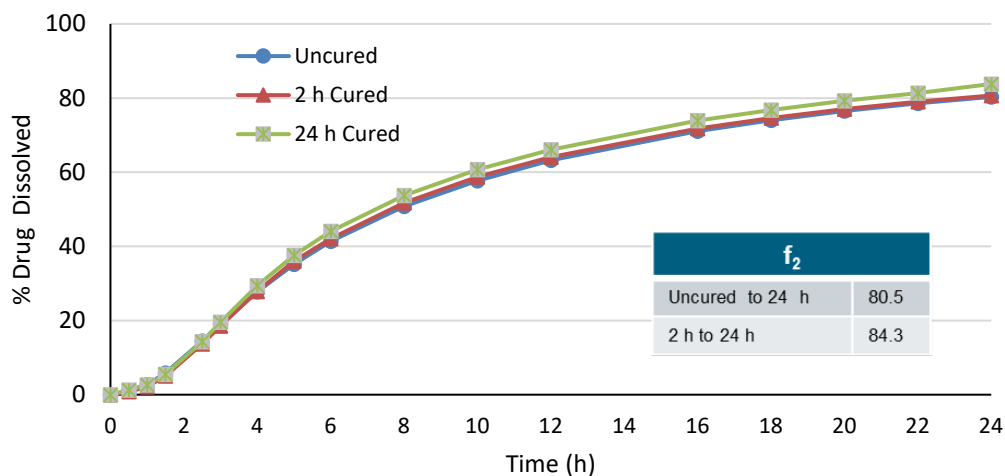



Figure 2. Drug Release from Propranolol HCl MR MPs Coated with Partially Formulated Ethylcellulose Dispersion (Aquacoat); Effect of Curing Duration

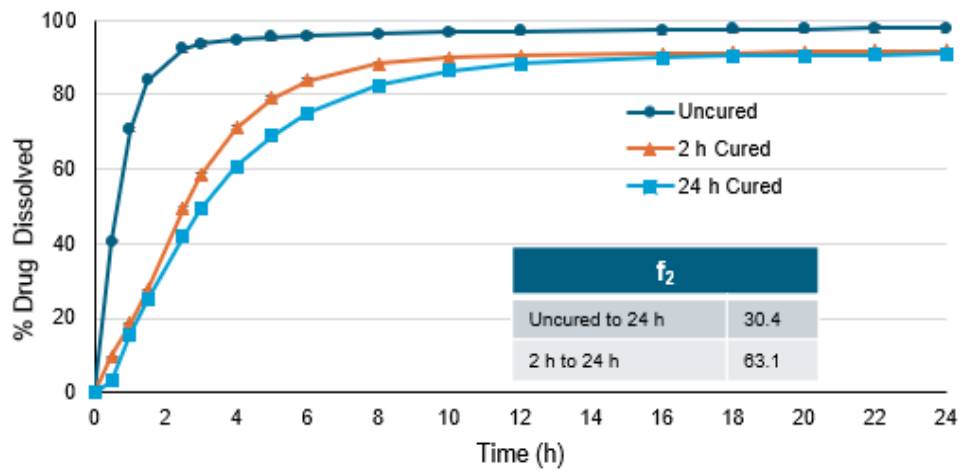

Results: SEM Images

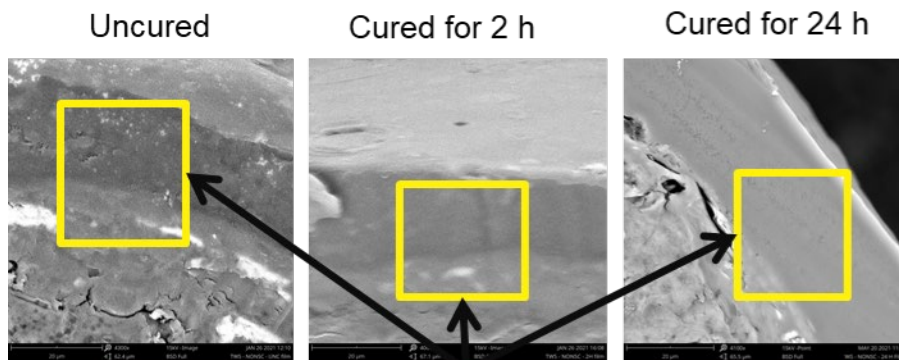

Surelease BM film remained dense even without curing

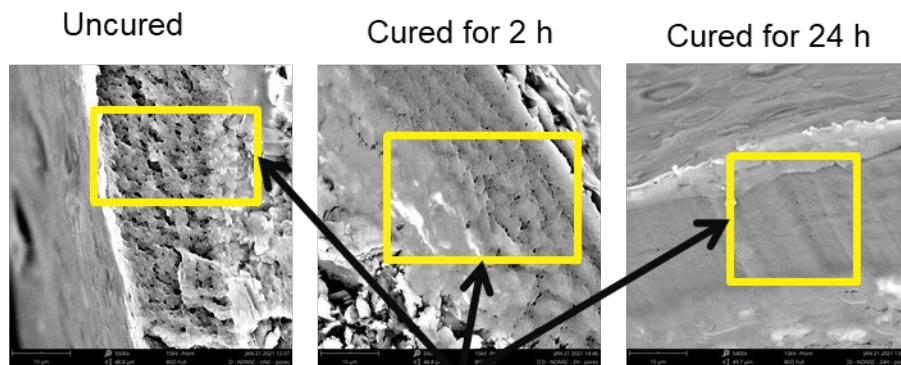

Aquacoat shows significant reduction in BM coating film porosity due to curing

Conclusions

Ethylcellulose aqueous dispersions were successfully coated to achieve controlled release for propranolol HCl multiparticulates. Consistent drug release was achieved for ethylcellulose aqueous dispersions. The results showed the importance of curing and processing parameters required for developing controlled release MP systems.

References

1. Moore, J. W. and H. H. Flanner, 1996, "Mathematical Comparison of Dissolution Profiles," Pharmaceutical Technology, 20 (6):64-74.

The information contained herein, to the best of Colorcon, Inc.'s knowledge is true and accurate. Any recommendations or suggestions of Colorcon, Inc. with regard to the products provided by Colorcon, Inc. are made without warranty, either implied or expressed, because of the variations in methods, conditions and equipment which may be used in commercially processing the products, and no such warranties are made for the suitability of the products for any applications that you may have disclosed. Colorcon, Inc. shall not be liable for loss of profit or for incidental, special or consequential loss or damages.

Colorcon, Inc. makes no warranty, either expressed or implied, that the use of the products provided by Colorcon, Inc., will not infringe any trademark, trade name, copyright, patent or other rights held by any third person or entity when used in the customer's application.

For more information, contact your Colorcon representative or call:

North America	Europe/Middle East/Africa	Latin America	India	China
+1-215-699-7733	+44-(0)-1322-293000	+54-1-5556-7700	+91-832-6727373	+86-21-61982300

You can also visit our website at www.colorcon.com



© BPSI Holdings LLC, 2021.

All trademarks, except where noted, are property of BPSI Holdings, LLC.

Aquacoat® is a trademark of International Flavors and Fragrances Inc. or its affiliates. © 2021 IFF. All rights reserved

CRS_2021_Mehta_Surelease