# The Effect of Film Coating on the Stability of Extended Release Metformin Hydrochloride POLYOX™ Matrices

Hua Deng, Scott Vass, Sandip Tiwari, Thomas P. Farrell and Ali Rajabi-Siahboomi

## **PURPOSE**

Polyethylene oxide (PEO) is prone to oxidation and thermal degradation at high temperatures, resulting in reduction of polymer viscosity and polymer chain cleavage during storage. PEO degradation can be regulated by oxygen impediment, oxygen removal and/or other methods as mentioned in the literature. The aim of this study was to investigate the effect of an Opadry If film coating system on the stability of metformin HCI extended released (ER) polyethylene oxide (POLYOX M) matrix tablets at accelerated storage condition (40°C/75% RH). Tablets were pulled at pre-determined time points and evaluated for physical properties and drug release. The effects of filler and presence of desiccant were also evaluated in this study.

## **METHODS**

#### **Materials and Formulations**

The composition of metformin HCl formulations is shown in **Table 1.4** Metformin HCl was used as the model drug, POLYOX Coagulant as the release controlling agent, and StarCap 1500® and/or microcrystalline cellulose (MCC) utilized as a filler in this study.

Table 1. Composition of Metformin HCI ER Formulations

Ingredients	% Composition (w/w)		
	F1	F2	F3
Metformin HCI (Wanbury Limited, India)	50.0	50.0	50.0
POLYOX Coagulant (IFF., USA)	20.0	20.0	20.0
StarCap 1500 (Colorcon Inc., USA)	29.0	-	14.5
Microcrystalline cellulose (Microcel 102SP, Blanver, Brazil)	-	29.0	14.5
Fumed silica (Aerosil 200, Evonik Industries, Germany)	0.5	0.5	0.5
Magnesium stearate (Akros Chemicals, Holland)	0.5	0.5	0.5
Total	100.0	100.0	100.0

## **Tablet Preparation and Film Coating**

The API, POLYOX Coagulant, filler(s) and fumed silica were passed through an ASTM 30 mesh (600  $\mu$ m), and then mixed in a 4 quart V blender (Patterson-Kelley Co., USA) at 25 rpm for 10 mins. Magnesium stearate was screened through an ASTM 40 mesh (400  $\mu$ m) then added to the powder mixture, followed by blending for an additional 3 mins.

Matrix tablets (1000 mg) were compressed on a rotary tablet press (Piccola, RIVA, Argentina), using caplet shaped tooling (19.1 mm  $\times$  9.3 mm), at the compression force of 30 kN (compression pressure of 190 MPa). <sup>4</sup> Tablets were then coated with a PVA-based Opadry II film coating system to the theoretical weight



gain of 4%, in a fully perforated coating pan (Labcoat I, O'Hara, Canada) fitted with a 10" pan without baffles. The film coating process parameters are shown in **Table 2**. Coated and uncoated matrix tablets were then packed in HDPE bottles, with or without a desiccant, and stored in 40°C/75% RH stability chamber for 6 months.

Table 2. Coating Process Parameters

Process parameters	Values	
Nozzle size (mm)	1.2	
Pan charge (g)	600	
Pan speed (rpm)	20	
Inlet air temperature (°C)	65-69	
Exhaust air temperature (°C)	48-53	
Product temperature (°C)	42-45	
Air volume [cfm / (m <sup>3</sup> /hr)]	130 / 218	
Atomization pressure (psi / bar)	15 / 1.0	
Pattern air pressure (psi / bar)	20 / 1.4	
Spray rate (g/min)	12-13	
Coating time (min)	10	

## **Tablet Characterization and Drug Release**

At each pre-determined time interval: tablet weight, breaking force, diameter and thickness were measured with an automated Multicheck tablet tester (Erweka, Germany). Tablet friability was also measured using a VanKel friabilator at 100 revolutions, 25 rpm (Varian Inc., USA). Drug release testing was performed using USP Apparatus II (VK 7000, Varian, USA) with sinkers at 100 rpm in 1000 mL of deionized water at 37 ± 0.5°C. Drug release was detected at a wavelength of 233 nm using a UV visible spectrophotometer (Agilent 8453, Agilent Technologies, USA) fitted with quartz flow cells of 1.0 mm path length.

### **Data Analysis and Model Fitting**

The release exponent (n) and release rate constant (k) were calculated by fitting the dissolution data to the Power Lawequation:

$$Q = k x t^n$$
 Equation (1)

where Q is the fractional amount released at time t, k is the kinetic constant, and n is the release exponent. In addition, the similarity factor ( $f_2$ ) was calculated by comparing the dissolution curves obtained at a given stability interval versus time zero results.<sup>5</sup>

### **RESULTS**

## **Physical Properties of Matrices**

The physical properties of the ER metformin matrices at time zero are shown in **Table 3**. All tablets exhibited acceptable mechanical strength, low weight variation and low friability values.



POLYOX™ -2-


The tensile strength of metformin ER matrices is shown in **Figure 1**. At time zero, the mechanical strength of the tablets was in the range of 0.86-1.34 MPa (uncoated) or 1.26-1.59 MPa (coated). At 6 month time point, uncoated tablets exhibited loss of mechanical strength (0.40-0.98 MPa), while coated tablets had minimal change (1.32-1.73 MPa) when packed with a desiccant. The enhanced mechanical strength and stability of coated tablets was attributed to the Opadry II film coating and its excellent moisture and/or oxygen barrier characteristics. The mechanical strength loss of uncoated tablets might be explained by the elastic nature of POLYOX and small moisture uptake. No significant differences in tablet weight or thickness were observed for uncoated or coated PEO tablets.

Table 3. Physical Properties of Metformin HCl PEO Matrix Tablets at Time Zero (Compression force/ Pressure of 30 kN/ 190 MPa) (n = 20)

Formulation	Opadry II	Hardness	Tensile	Tablet	Weight	Friability*
(Filler)	Coating	(kp)	Strength	Thickness	Variation	(%)
			(MPa)	(mm)	(%)	
F1 (StarCap	No	10.40 ± 0.75	$0.86 \pm 0.07$	6.76 ± 0.02	0.44	0.49
1500)						
F2 (MCC)	No	15.70 ± 0.56	1.34 ± 0.05	$6.56 \pm 0.02$	0.36	0.30
F3 (StarCap	No	14.30 ± 0.53	1.24 ± 0.05	$6.47 \pm 0.02$	0.48	0.49
1500+MCC)						
F1 (Starcap	Yes	15.90 ± 1.30	1.26 ± 0.10	$7.03 \pm 0.04$	0.71	0.00
1500)						
F2 (MCC)	Yes	19.50 ± 1.00	1.59 ± 0.08	6.85 ± 0.04	0.70	0.00
F3 (StarCap	Yes	18.30 ± 1.10	1.51 ± 0.09	6.94 ± 0.02	0.46	0.00
1500+MCC)						

Note: \*n=10 for tablet friability testing.

Figure 1. Mechanical Strength of Metformin HCI ER Matrix Tablets (n=20): (A) No Desiccant (B) With Desiccant

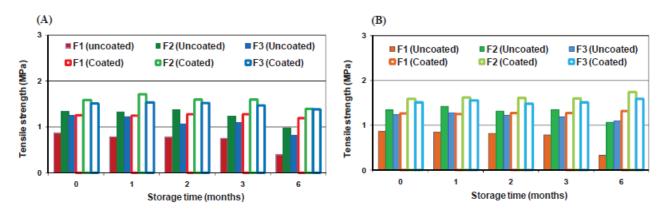

## **Drug Release Profiles**

Metformin HCI release profiles for uncoated and coated tablets are shown in **Figures 2-4**, respectively. Similar drug release profiles ( $f_2 > 50$ ) were obtained for tablets at time zero and at 6 months regardless of filler choice or film coating. The Opadry II film coating system had no impact on drug release profiles ( $f_2 > 65$ ). All formulations showed good data fitting to Power Law equation ( $R^2 > 0.99$ ), and the release exponent (n) was in the range of 0.59-0.61 indicating drug release was mainly controlled by diffusion.<sup>6</sup>

Figure 2. Metformin HCl Release Profiles from PEO Tablets with StarCap 1500 as Filler (n = 6): (A) Uncoated, No Desiccant ( $f_2 = 65-85$ ); (B) Coated, No Desiccant ( $f_2 = 72-73$ )

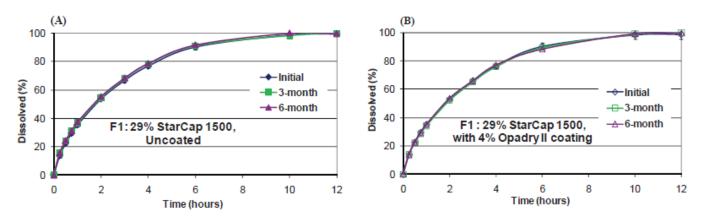
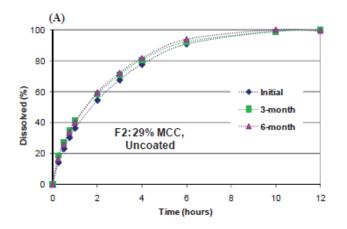


Figure 3. Metformin HCl Release Profiles from PEO Tablets with MCC as Filler (n = 6): (A) Uncoated, No Desiccant ( $f_2 = 67-78$ ); (B) Coated, No Desiccant ( $f_2 = 53-79$ )





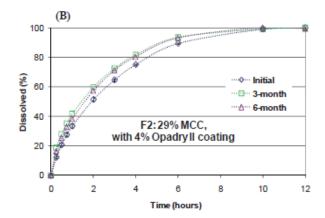


Figure 4. Metformin HCl Release Profiles from PEO Tablets with StarCap 1500 + MCC as Filler (n = 6): (A) Uncoated, No Desiccant ( $f_2 = 79-82$ ); (B) Coated, No Desiccant ( $f_2 = 78-83$ )

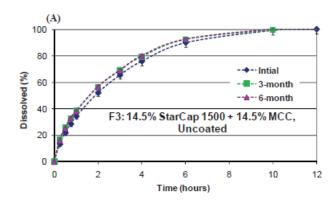

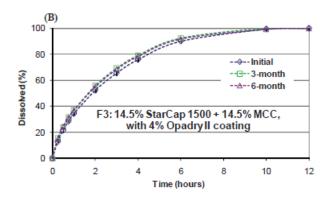

### CONCLUSIONS

The application of PVA-based Opadry II film coating system, at 4% weight gain on metformin HCI PEO matrices, significantly enhanced tablet stability and prevented loss of tablet mechanical strength during accelerated storage conditions. This is due to the excellent moisture/oxygen barrier of PVA based coating systems. The presence of a desiccant in the package provided additional protection and led to consistent mechanical strength values during storage. Drug release from the matrix tablets was not significantly affected by the presence of Opadry II film coating or various fillers. Similar drug release profiles were obtained for all evaluated matrices upon storage at an accelerated condition for 6 months.

## **REFERENCES**

- 1. Voorhees K.J., et al., An investigation of the thermal degradation of PEG, *J. Anal. Appl.* Pyrol. 1994;30 (1): 47-57.
- 2. Rajabi-Siahboomi A.R., & Farrell T.P., The applications of formulated systems for the aqueous film coating of pharmaceutical oral solid dosage forms, Editors: Linda Felton and James McGinity, *Aqueous polymeric coatings for pharmaceutical dosage forms*, 3rd Edition (2008).
- 3. Palmer, D., et al., The Influence of Film Coating and Storage on Propranolol Hydrochloride Release from Polyethylene Oxide Extended Released Matrices. 38th Annual Meeting & Exposition of the Controlled Release Society, National Harbor, MD, 2011.
- 4. Deng, H., et al., Evaluation of StarCap 1500® as an Excipient in Extended Release Hydrophilic Matrix Systems. AAPS annual meeting and exposition, New Orleans, LA, 2010.
- 5. Moore J.W., & Flanner H.H., Mathematical comparison of dissolution profile, *Pharm. Tech.* 1996;20: 64-74.
- 6. Korsmeyer, R.W., et al., Mechanism of solute release from porous hydrophilic polymers. *Int J Pharm* 1983;15: 25-35.

The information contained herein, to the best of Colorcon, Inc.'s knowledge is true and accurate. Any recommendations or suggestions of Colorcon, Inc. with regard to the products provided by Colorcon, Inc. are made without warranty, either implied or expressed, because of the variations in methods, conditions and equipment which may be used in commercially processing the products, and no such warranties are made for the suitability of the products for any applications that you may have disclosed. Colorcon, Inc. shall not be liable for loss of profit or for incidental, special or consequential loss or damages.

Colorcon, Inc. makes no warranty, either expressed or implied, that the use of the products provided by Colorcon, Inc., will not infringe any trademark, trade name, copyright, patent or other rights held by any third person or entity when used in the customer's application.

For more information, contact your Colorcon representative or call:

North America Europe/Middle East/Africa +1-215-699-7733 +44 (0)-1322-293000

Latin America India China +54-11-5556-7700 +91-832-6727373 +86-21-61982300



© BPSI Holdings LLC, 2016

The information contained in this document is proprietary to Colorcon, Inc. and may not be used or disseminated inappropriately.

POLYOX™ is a trademark of International Flavors and Fragrances Inc. or its affiliates. © 2021 IFF. All rights reserved

All trademarks, except where noted, are property of BPSI Holdings LLC.